


SMART FITNESS SCALE

Quick Start Guide



Install "Starfit" and complete registration with email address or login your existing account.

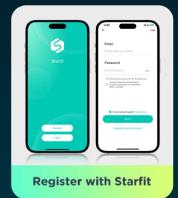

Ensure to keep Bluetooth ON; and for Android system, GPS & Location must be ON along with Bluetooth.



Find "My Devices" in the Starfit APP, click "+" in the upper right corner to search and connect the scale (Bluetooth name: Runstar-Rx).



Open the "Measure" screen on the App, step on the scale to start measuring. The App will automatically read and record all results.

Trouble Shooting

Why my App cannot connect with my scale?

- 1. Make sure the network on your phone works well.
- 2. Make sure the Bluetooth is on and GPS & location is allowed for Android system.
- 3. Make sure the scale is on while pairing.
- 4. If you cannot connect with the App, please follow below steps.
- 5. You may have closed some request authorization. Please reinstall and accept all authorizations and log in to verify again.



Click "Account"



Click "Device"



Click "+" to search the Bluetooth device



Return to Homepage, it shows "connected" when paired

Why I cannot check my body fat rate?



No shoes, no socks on. Measure with bare feet.



Make sure the bare feet touch the metal sensors.



Stand on the scale for at least 5 seconds.



Make sure your App is on the homepage when measuring.

Why my scale is inaccurate?





Customer Support

service@runstar.store