# Digital Weight Scale

## USER MANUAL

**Model: BG2403-B**

# Contents

---

|           |       |
|-----------|-------|
| • ENGLISH | 01-03 |
|-----------|-------|

---

|           |       |
|-----------|-------|
| • DEUTSCH | 04-06 |
|-----------|-------|

---

|            |       |
|------------|-------|
| • FRANÇAIS | 07-09 |
|------------|-------|

---

|            |       |
|------------|-------|
| • ITALIANO | 10-12 |
|------------|-------|

---

|           |       |
|-----------|-------|
| • ESPAÑOL | 13-15 |
|-----------|-------|


Thank you for purchasing the RunStar Digital Weight Scale.  
To ensure safe and proper use of this product, please read all instructions and safety guidelines in this manual carefully.  
Keep this manual for future reference.

## Specifications


1. Product Name: Digital Weight Scale
2. Model: BG2403-B
3. Weight Capacity: 400lb / 180kg
4. Increments: 0.1lb / 50g
5. Weight Units: lb / kg
6. Display: LED
7. Battery: 3×1.5V AAA Batteries (inserted)
8. Dimensions: 280×280×20 mm
9. Platform: Tempered Glass
10. Auto-off: 10 Seconds

## Instructions for use

1. Remove the plastic strip from the battery compartment before first use.

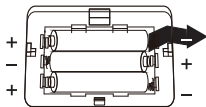2. Place the scale on a hard, flat and level floor.  
Avoid using the scale on carpets or soft surfaces.

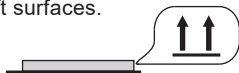
3. Step on the scale and then step off to activate the device.  
Wait for the screen to display "0.0".
4. Step on the scale and stand still until the display flashes, showing the final measurement.
5. The scale will automatically power off after 10 seconds of inactivity.
6. To switch units (kg-lb): Press the "UNIT" button on the back of the scale display.

## Notice

1. For consistent results, it is recommended to take measurements at the same time each day.
2. Replace the batteries when the low battery sign “LO” appears on the LED display.

3. If “Err” appears on the display, it indicates the scale has been overloaded.
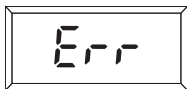
## Caution

1. Remove your shoes before stepping on the scale to avoid scratching or damaging the tempered glass platform.
2. Ensure that both the scale platform and your feet are dry before use. The glass platform is slippery when wet.
3. Do not jump or stomp on the scale.
4. Clean the scale with a soft, dry cloth or a slightly damp cloth if needed. Avoid using chemical or abrasive cleaners, as they may damage the surface.
5. Keep the scale away from water and heat sources to prevent damage.
6. The scale is intended for home use only.

## Package Contents

- 1 × RunStar Digital Weight Scale
- 3 × 1.5V AAA Batteries
- 1 × User Manual

### ***FCC Statement***

Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference.
- (2) This device must accept any interference received, including interference that may cause undesired operation.

**FCC Statement NOTE:** This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

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