SMART BODY FAT SCALE FG2019LB

User Manual

Contents

• ENGLISH	01-09
-----------	-------

• DEUTSCH	10-18
-----------	-------

• FRANÇAIS	19-27
------------	-------

• ITALIANO	28-36
------------	-------

• ESPAÑOL	37-45
-----------	-------

Bluetooth Scale Body Index

- Body weight
- Body balance tracking
- BMI (Body mass index)
- Body fat rate
- Fat-free body weight
- Subcutaneous fat
- Heart rate
- Visceral fat
- Body water
- Skeletal muscle
- Muscle mass
- Bone mass
- Protein
- BMR
- Body age
- Cardiac Index

Trouble shooting

1. Battery volume
2. Connected
3. Weight units
4. Body fat rate
5. Body type assessment
6. Body weight
7. BMI
8. Body heart rate

No display

Change battery

Change battery

Overload

Calibrating

For body weight: Place the scale on an even surface and wait at least 10 seconds before stepping on the scale. The reading should display 0.0 kg before stepping on.

During body fat testing process, the display shows bars in different colors in sequence, which refer to different body types.
After measured, the display shows the body fat rate and determinate color bar.

The correspondence between colors and body types:

deep blue=thin


light blue=normal-

green=normal+

orange=fat

red=obese

For example: The below display shows the **body ⑤ types** is normal.

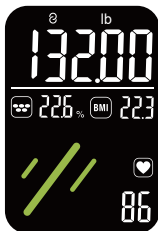Caution

Conductive Areas

1. Do not shake or drop the scale.
2. Do not immerse the scale in water or use chemical cleaning agents. Clean the scale with a slightly damp cloth.

3. Do not place the scale on any carpet or uneven flooring. Always place the scale on hard, flat, dry surface before measurement.

4. Always keep your scale in horizontal position.

5. Any measurement obtained using this device is for reference only and should not be considered as a medical opinion.

6. Consult your primary care doctor or physician before making changes to your diet, exercise plan or physical activity.

Specifications

Size: 11.02*11.02*0.98in/280*280*25mm

Weight: 3.44lb / 1.56kg

Batteries: 4 x 1.5V AAA batteries

Capacity: 400 lb / 180 kg / 28 st

Division: 0.1lb/50g

Unit: lb/st/kg

Accuracy Range

50kg:±0.3kg	100kg:±0.4kg	180kg:±0.5kg
110lb:±0.66lb	220lb:±0.88lb	400lb:±1.1lb


Operating instructions

The Bluetooth Scale can connect with your mobile phone (Android 6.0+ & IOS 8.0+) or iPad via Bluetooth 4.0. Download the 'Starfit' app free from the App Store or Google Play.


Pairing the scale with your smartphone


Turn on Bluetooth on your smartphone.

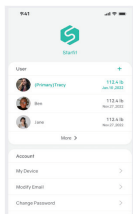Press [Device] and [Device Management] to enter measuring page and find device.

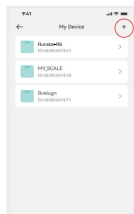
Power on the scale, press[+] to select "Runstar-R6" to add the device.



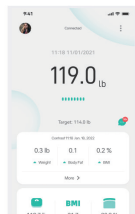
Come back to measure page. Step on the scale and keep standing for 5-10 seconds with bare feet. The weighing has finished after the analysing data come out.

View reports, delete data, set reference data, share progress by recently, weekly, month and year.

Click [Chart] to view your progress history. Enter it into the calendar to form a list and select a dataset or all data you want to delete. You can also select different date data to contrast and share to friend via facebook, instagram, or SMS .



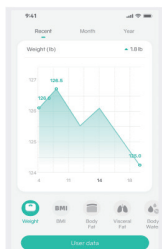
At the bottom, you will find a bar where you can select different types of data to display.

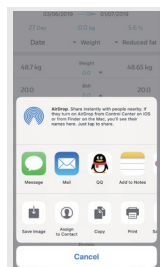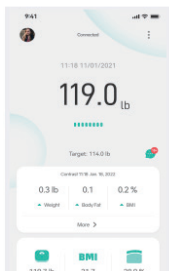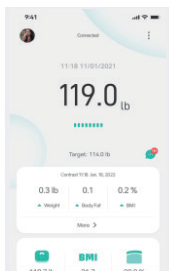Add/Delete Users & Basic Settings

Under [Account] you can add user(s) [+] . To delete an account, swipe the account to the left and select [Delete Account] . Themes, Set goals, Weight units, Weighing done sounds, Languages, Passwords can all be set under the Account page. Please note that you can change the weighing unit by pressing the Unit button. , or you can also switch units from the Setting page of App.

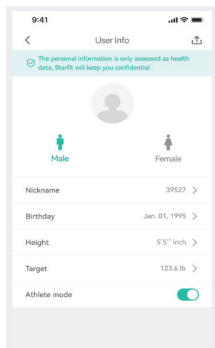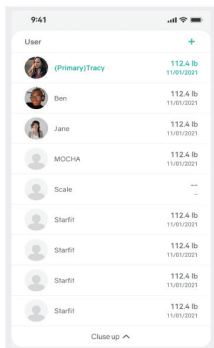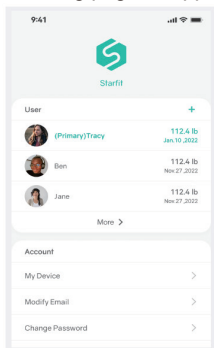Frequently Asked Questions&Answers

1.Why doesn't the scale work? Why do the data on screen disappear in a flash?

- Please check if the batteries are installed properly, replace batteries if necessary.

2. The scale cannot connect with App.

- Check if the phone's software is iOS 8.0 or higher or Android 6.0 or higher.
- Download and open the newest version of the Starfit App.
- Enable Bluetooth on your phone. For Android 6.0 or higher, location needs to be enabled.

3.No body fat measured when weighing.

- Step on with dry and bare feet.
- Ensure Bluetooth is enabled and working.
- Ensure personal data is entered.

Warning Note

Changes or modifications not expressly approved by the party responsible for compliance could avoid the user's authority to operate the equipment.

Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference.
- (2) This device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

CUSTOMER SERVICE

WhatsApp (chat only): +1 (626) 743-7007

Email: service@runstar.store

Monday - Friday

FCC ID: 2AP3Q-FI2019LB-B

Model No.: FG2019LB

Guangdong Welland Technology Co., Ltd.

No.85, Minke East Road, Mining Science & Technology Park,
Shiqi District, Zhongshan, Guangdong, China

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